



E Light

YB 8000

The latest in RF and IPL combination technology.

1 Year warranty on all equipment
Extended warranties available
Nationwide servicing
Training provided

Treatment Principle

E-light machine combines IPL and RF so it can improve the client treatment experience through a significant reduction in pain levels. It combines three advanced technologies:

- Bipolar Radio Frequency
- Intense Pulse Light
- Skin Contact Cooling

When the three are united in one treatment client results are significantly improved. The energy of RF can reach middle layer of the skin and heat tissue, thus lower energy is applied during IPL treatment. Pain felt during IPL treatment will be significantly reduced and far better results will be achieved in all applications. In addition the cooling system assists in further client contentment. The absorption of the RF energy is not concerned with melanin so the E-light treatment can obtain excellent results on deep colour skin (skin treatments) and light hair (hair removal treatment). Special designed ISL device can monitor the skin temperature to insure safety of treatment and smart technology of hand piece allows direct contact of tip onto skin.

Handles

One E-light Handle 12mmX45mm with 7 filters : 420nm, 520nm, 560nm, 585nm, 640nm, 690nm, 755nm

Treatments

- skin tightening: skin lifting, improve flexibility and tone of skin
- skin rejuvenation: fine wrinkles, pigments, freckles, vascular lesions (telangiectasis) removal, skin whitening, shrink large hair pores
- hair removal: all skin types and hair colour
- acne removal: improve the condition of oily skin, kill acne bacil

One RF mono polar handle with 5 heads (preferred for optimal body skin treatment and fat reduction)

One extra RF bipolar handle with 2 heads (preferred for optimal facial skin treatment)

General functions of RF

- deep wrinkle
- skin lifting
- promote metabolism of cell
- correct facial contour
- improve skin elasticity
- improve ATP, treatment of acne
- diminish inflammation



PREFMED

PREFERRED MEDICAL TECHNOLOGIES

Distributed by Pref Med
Tel: +27 021 791 2870 Fax: +27 021 791 2871
www.prefmedtech.com



ISO Certification No. 07 0142 SJ
Department of Health License No.790/15175



E Light


YB 8000

The latest in RF and IPL combination technology.

Specifications

Wavelength	Changeable--420/510/560/585/640/690/755-1200nm
IPL Energy Density	10-50J/cm2
Pulse Sequences	1-17 pulses adjustable
Pulse Delay	1-50ms
Pulse Duration	1-20ms
Repeat Frequency	0.3-1Hz
RF Frequency	1, 5, 10 MHz
RF Energy Density	1-100J/cm3
Light Spot	12mmX45mm
Display	8.4 inch colour touch screen Bigger display (8.4 inch), , more convenient for the
IPL Power	1200W
RF Power	1000w
Power Supply	110V/220V

Before And After
E-light / IPL

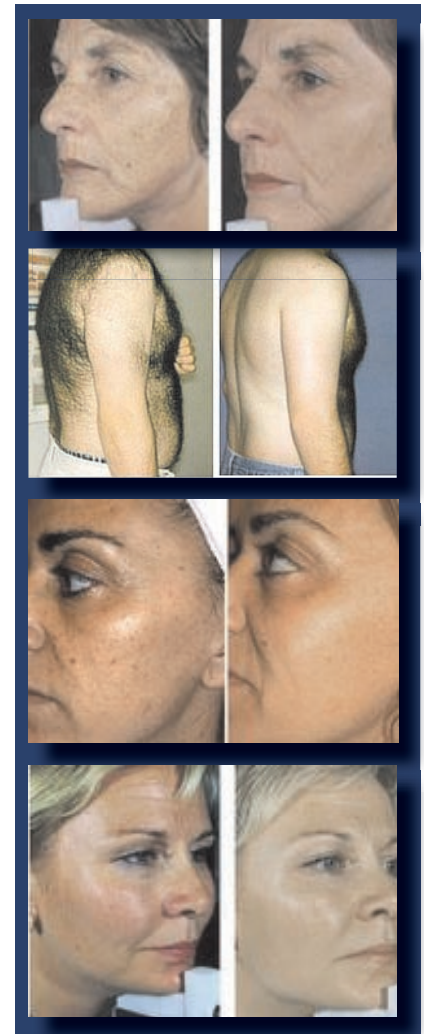

Before And After
E-light / RF



PREFMED

PREFERRED MEDICAL TECHNOLOGIES

Distributed by Pref Med
Tel: +27 021 791 2870 Fax: +27 021 791 2871
www.prefmedtech.com



ISO Certification No. 07 0142 SJ
Department of Health License No.790/15175